

THE KIT CONTAINS:

Tirol Kit cod. 21-0120-023

In just 10,8 Kg

1 LIGHT MULE
ultra-light single-wheel system

4 HANDLES
handles for Titan T

2 SS-SHOULDER SUPPORT
lifting harnesses

1 BAG TIROL KIT
storage bag

Handles can be adjusted in **12 DIFFERENT POSITIONS** as needed, with short or long versions enabling the stretcher to fit through even the narrowest of places in complete safety.

The handles can be positioned to function as legs for the frame, providing a stationary position for the stretcher and enabling the casualty to be attended to whilst **COMPLETELY LIFTED OFF THE GROUND.**

ShoulderSupport: shoulder harnesses which are easy to fit, adjustable and are fitted with a high visibility system, enabling more efficient distribution of the weight to be transported.

LOW WEIGH (total weight of kit 10.8 kg), which can be distributed among several operators and is backpackable, taking up minimal space.

PATENTED LOCKING SYSTEM

compliant with regulation 93/42

TITAN PROTECTION.

Special coating which protects the Titan-T basket stretcher enabling it to be safely dragged over all kinds of terrain.

The whole kit can be put together in a matter of seconds by a **single operator** and the transport wheel can be fitted whilst the casualty is in position.

DIVISIBLE & BACKPACKABLE

Load capacity 225 Kg
the lightest on the market only 7 Kg

Light Mule: single-wheel system comprising a foldable "V" support and fast clip lock. Equipped with a high load capacity wheel with tyre designed for mountain terrain, offering high adherence and manoeuvrability off-road as well.


TITANIUM FRAME

Anti-wear and ant-abrasion treatment of metal parts with Niploy, ensuring durability over time.



basket stretcher Titan-TI

- ▶ Extremely light- due to its 100% titanium frame and TIG welds, its weight is reduced by half compared to stainless steel stretchers.
- ▶ The high-density polyethylene netting prevents rust from forming, ensuring greater sturdiness and durability over time.
- ▶ The top tube (25.4 mm) gives a safer and more comfortable grip, especially in the event of long journeys.
- ▶ With its rapid closure system you can put the stretcher together and take it apart in just 15 seconds.
- ▶ Using "StratLoad™" attachment points the stretcher can be quickly connected to a litter bridle for lifting. These attachment points ensure safe transport without risk of damaging the frame.

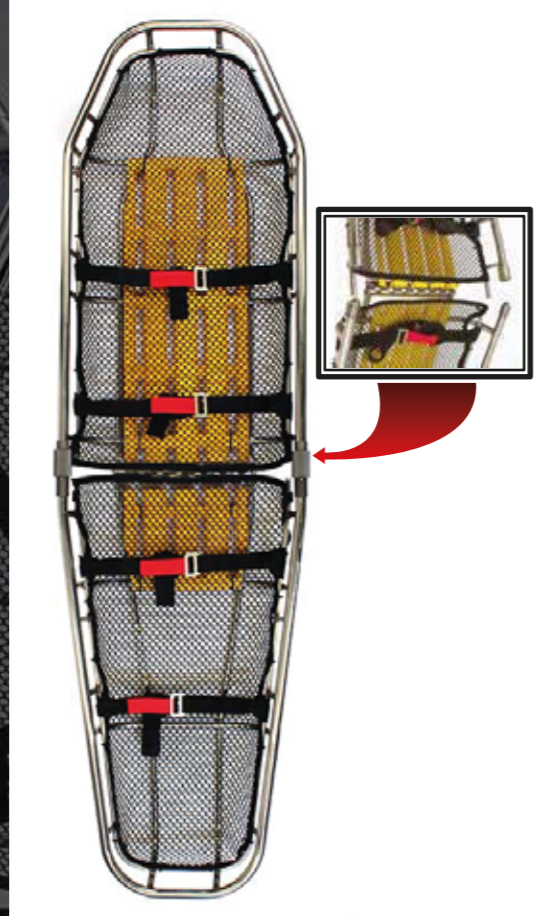

CERTIFICATIONS:
 ▶ CE CLASSIFIED
 ▶ CERTIFIED CLASS 1



CODE	DESCRIPTION	WIDTH	LENGTH	DEPTH	LOAD	COEFFICIENT
11-0119TIROL	Titan TI	60 cm	212 cm	18 cm	7 kg	1134 kg

developed in
collaboration with
**Tirol
Bergrettung**

OPTIONAL ACCESSORIES

Titan Protection
cod. 21-0120-010
Titan Protection (divisibile)
cod. 21-0120-010_1



For further information
visit www.ferno.it

Ears (sailing Anti-rotation)
cod. EARS



For further information
visit www.ferno.it

Litter bridles
cod. 15-0108 (hard)
cod. 15-0109 (soft)



For further information
visit www.ferno.it

Tirol kit bag
cod. BAG TIROL KIT



For further information
visit www.ferno.it

Travers rescue
(travers rescue)
cod. 74000



For further information
visit www.ferno.it

consult us with
QR CODE READER



Numero Verde
800 501711

FERNO

Reg. Office & Sales
Via B. Zallone n. 26
40066 Pieve di Cento [Bo] | Italy
P +39 051 6860028 | F +39 051 6861508
info@ferno.it
www.ferno.it

Just imagine how far
you could get with
only 11 kilograms.

FERNO