

OPTIONAL ACCESSORIES

MonoClack
Quick locking system part
cod. 65-0430-001

MONOCLACK
TESTED **10G**



For more information visit
www.ferno.it
www.ferno-schweiz.ch

Battery
cod. MWK-BATTERYM28



For more information visit
www.ferno.it
www.ferno-schweiz.ch

FERNO

Click, Touch & Go!
Monorudder

consult us with
QR CODE READER



FERNO
Italia

Sede Legale e Uffici Commerciali
Via B. Zallone n. 26
40066 Pieve di Cento [Bo] | Italy
P +39 051 6860028 | F +39 051 6861508
info@ferno.it
www.shop.ferno.it  **800 501711**

FERNO
Schweiz

FERNO (Schweiz) AG - Cheerstrasse 16,
6014 Littau/Luzern - Switzerland
P +41 41 259 60 00
F +41 41 259 60 01
MWST-Nr.: CHE-100.510.974
www.ferno-schweiz.ch



MONORUDDER

MOTORIZED INTERFACE FOR NON-RECLOSABLE TRANSFER CHAIRS



MOTORIZED

The powerful motor allows the transport of patients on all types of stairs with ease, without effort and in total safety.

TECHNICAL SPECIFICATIONS

WEIGHT	22 KG	CAPACITY	230 KG
HEIGHT	107 CM	POWER SUPPLY	28V LI-ION BATTERY
WIDTH	42 CM	MATERIALS	ALUMINIUM
DEPTH	22 CM	IP RATING	IP56

QUICK AND INTERCHANGEABLE



1 Lift the Monorudder.



2 Insert the pin into the proper grooves.



3 Press the device towards the chair backrest.

COMPATIBLE WITH CHAIRS


Dlouhy Utila Venice pick up



SLIM

Thanks to its small size, Monorudder can be stored in very little space.

The additional rear wheels make it portable on its own even when not attached to the transfer chair.

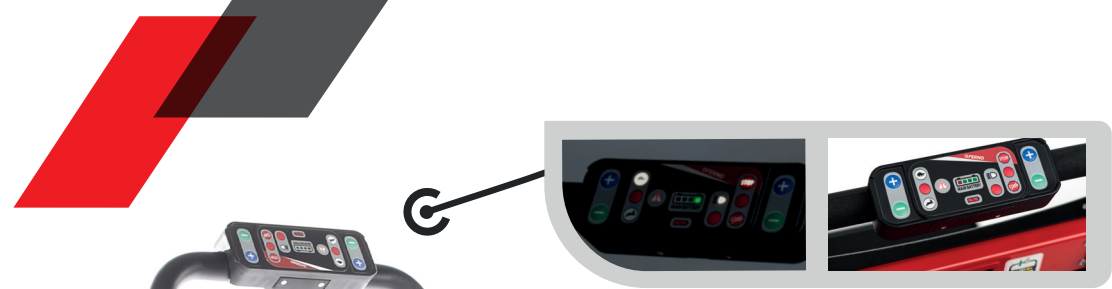

REMOTE CONTROL



The new control system guarantees use even in the most extreme temperature conditions. Shockproof and Water Resistant. Resistant to dust and water infiltrations (IP56 certified).

Backlit controls



ANTIURTO



COMANDI TOUCH PROGRESSIVI



COMANDI RETROILLUMINATI



WATER RESISTANT

ADJUSTABLE HANDLES AND ADJUSTABLE HEIGHT



SAFETY LIGHTING

To provide visibility in dark conditions (poorly lit stairways or with a timed light), thus ensuring transport in total safety.

