



Revolutionary device for extrication and spinal immobilisation. Designed to face up to difficulties in restricted environments, protecting the vertebral column of the patient. Its way of use significantly reduces the application time and at the same time increases effectiveness.



**SAFETY**

Carbon stiffness and lightness allow a quick and safe extrication also in confined spaces and under harsh conditions.



**EASE**


The coloured buckles facilitate and speed up the patient immobilisation process.

**EXTRICATION IN LESS THAN 20 SECONDS**

**XT = CARBON BOARD + CLASSIC RESTRAINT SET + QHI + CARRY BAG + WIZ LOC COLLARS (3PCS)**



**part no. Kit XT PLUS-B**

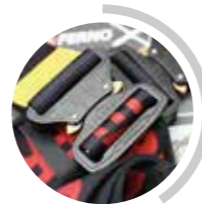

**COBRA BUCKLE**

Capacity 9 kN (917 kgf)  
Aluminium alloy 7075

**LEG HOOPS**



**LOOPER DY**

Capacity in **22Kn**



**Human LIFT**



▶ RESTRAINT SYSTEM TESTED IN ACCORDANCE WITH EN1498




**XT Pro = XT + QHI + Pro restraints + Carry bag (SEE BACK PHOTO)**

**part no. XT PRO**




Xt floating is a kit designed to ensure extreme manoeuvrability and reliability during water rescue. The head float and the two side floats ensure an optimal patient support, preventing capsizing and protecting him or her from accidental impacts.

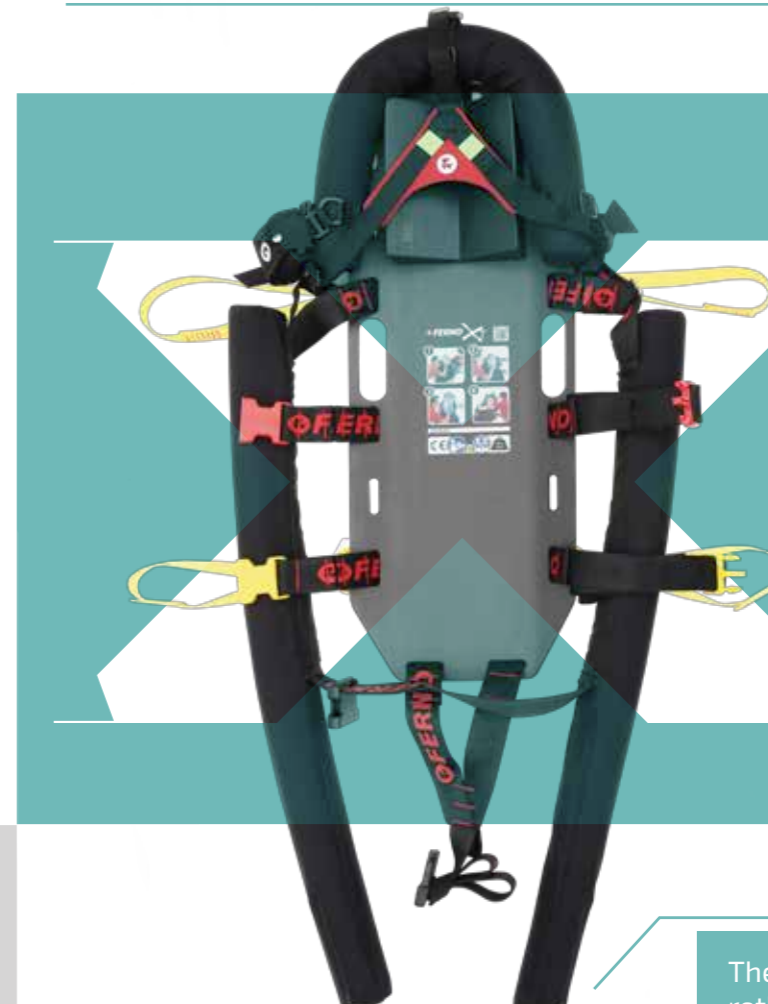

**FLOATING**

The layout of semi-rigid polyurethane floating tubes is conceived to **keep the patient with the head out of water**. Buoyancy is greater than the one of a buoyancy vest.

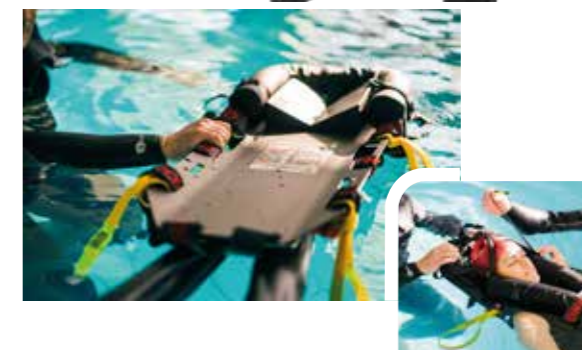

The HUMAN LIFT kit allows vertical rescue in isolated environments. The device keeps the patient posture on the rescue line, considerably reducing the handling space. Configuration speed is ensured by the single tape that facilitates patient lifting.

Quick float locking and unlocking thanks to **FAST CLIP** latches.

**FAST CLIP**



**XT Floating = XT + QHI + Floating restraints + Carry bag (SEE BACK PHOTO)**

**part no. XT Floating**

# Technical specifications



## XT

CARBON BOARD + Classic restraint set + QHI + Carry bag + Wiz Loc collars (3PCS)

Width	11.81 in - 300 mm
Length	32.68 in - 830 mm
Height	*2.36 in - 60 mm
Weight (including accessories)	4.40 Lbs - 2 Kg
Load capacity	352 Lbs - 160 Kg

\* height including restraints

## XT Pro

XT + KIT PRO

Width	11.81 in - 300 mm
Length	32.68 in - 830 mm
Height	*2.36 in - 60 mm
Weight (including accessories)	7.49 Lbs- 3.4 Kg
Load capacity	352 Lbs - 160 Kg

\* height including restraints

## XT Floating

XT + KIT FLOATING

Width	13.78 in - 350 mm
Length	47.24 in - 1200 mm
Height	*4.33 in - 110 mm
Weight (including accessories)	7.27 Lbs - 3.3 Kg
Load capacity	352 Lbs - 160 Kg

\* height including restraints



**DO YOU ALREADY OWN AN XT? UPDATE IT!**

*Our products are designed to be UPDATEABLE TO THE LATEST VERSIONS*

HERE ARE PART NUMBERS THAT MAY BE SUITABLE FOR YOU



### kit PRO

part no. KIT PRO

includes: 4 restraints with cobra buckle  
2 handling rings  
2 paddings for restraints (leg hoops)  
1 snap-hook  
1 lifting tape



### kit FLOATING

part no. KIT FLOATING

includes: 4 handling rings  
3 floating tubes  
2 transverse retaining strips  
1 longitudinal retaining strip



#### Italy

P +39 051 6860028  
F +39 051 6861508  
info@ferno.it

www.ferno.it



#### Schweiz

P +41 41 259 60 00  
info@ferno-schweiz.ch

www.ferno-schweiz.ch



Extrication  
for professionals

